

A-dec 2000



a-dec®

# The A-dec Difference.

The best way to discover what dentists need in a chair and delivery system is to ask. We did. And the A-dec 200 design reflects the way dental teams work, and purposefully supports efficiency in your practice. Everything about the A-dec 200 reflects the skill and quality of your care.



Backed by a 5-year warranty.





A dental unit featuring a white ceramic sink. Above the sink is a rectangular control panel. Below the sink, a white tray holds several dental instruments: a water spray nozzle, three handpieces (two silver and one grey), and a control console with a blue display and buttons. The brand name 'a'idec' is visible on the control console.

Designed for Durability.



**Pure light. Pure choice.** Two optional LED lights let you choose the best fit for your practice.

Your equipment needs to work when you do. That's why A-dec manages the manufacturing processes, so we can provide doctors with the most durable, reliable equipment. For example, we build and test each chair cylinder from raw material to the finished assembly, and back it with a 10-year warranty. And we manufacture our own high-quality delivery system tubing. Our production systems are consistently monitored to ensure that they meet A-dec quality standards. With the A-dec 200 solution, you get a strong, robust system designed to withstand years of constant use.





# Designed to Perform.

A-dec 200 is a system that's not only made to last, but made to grow and evolve with your practice. Detailed features such as a stain-resistant cuspidor that rotates toward the patient, a telescoping assistant's arm for close work, and a touchpad that integrates chair, light, and cuspidor functions all add up to a better workflow. For future expansion, an open platform leaves room to add or change ancillary devices. You get peak performance now, and years from now.



- Three axes light head lets you position the light at just the right angle for you and your patient.
- Four preset chair positions give you added flexibility for positioning the patient.
- Ancillary accommodation allows you to add new technology as your practice evolves.
- The chair's two-position armrests and seat articulation provide patient comfort.
- Sleek covers keep surfaces clean.

## Delivery System

Keep everything you need within reach to increase the speed of treatment. For added accessibility, position the adjustable instrument holders at the most comfortable angle.



A-dec 200



A-dec 300 Continental



A-dec 300 Traditional

## Touchpad

Fingertip access simplifies chair adjustments and provides quick access to dental light and cuspidor control. Optional deluxe touchpad adds controls for integrated clinical devices such as electric motors, cameras, and scalers.



Standard



Deluxe

## Foot Controls

Disc style with chip blower—or new A-dec lever style for precise modulation of electric and pneumatic handpieces.



Disc foot control with chip blower



Lever foot control

## Dental Lights

A-dec dental lights provide true light for accurate diagnosis, reduced eyestrain and unparalleled ergonomics.



A-dec 200 Halogen



A-dec 500 Halogen



A-dec 300 LED



A-dec 500 LED

## Monitor Mount

The monitor mount positions easily over the patient for convenient patient consultation and education.



*Monitor Mount*

## Cuspidor and Support Center

The porcelain cuspidor swivels for just the right patient position. Easy-to-fill, two-liter water bottle provides a continuous water supply through multiple procedures. The support center offers ample storage room to add ancillaries such as amalgam separators and vacuum systems.



*Cuspidor and support center*

## Assistant's Instrumentation

Three-position holder and telescoping arm conveniently positions vacuum instruments for patient access.




*A-dec (aluminum) with touchpad*



*Durr (polymer) without touchpad*

## Seamless Upholstery


A sleek vacuum-formed fit with no visible seams make cleaning and maintenance easy. See your authorized A-dec dealer for accurate color samples and the most current product information.

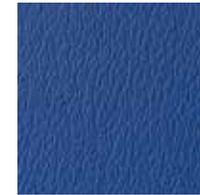


*sky blue*



*indigo*

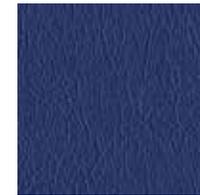

*sapphire*

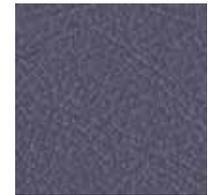

*amethyst*



*lemongrass*



*black*



*paprika*



*tawny*

## Specifications

### A-dec 200 Dental Chair

|                                  |  |
|----------------------------------|--|
| Range of lift height             | 15.75" to 28.75" (400 mm to 730 mm)  |
| Flexible backrest                | 2.25" (55 mm), seamless  |
| Drive system                     | Hydraulic base and tilt  |
| Chair back and base limit switch | Yes, stops downward movement when obstruction is detected. Patient comfort toeboard articulation |
| Headrest                         | Dual-articulating and gliding (knob style)   |
| Armrest                          | Pivoting   |
| Integrated 300-watt power supply | Yes  |
| Touchpad or foot switch          | Both optional (touchpad available on delivery systems)   |
| Mounting capability              | Floating - 2" post mount   |
| Upholstery options               | Seamless   |

### Delivery Systems

|                                     | <b>A-dec 200</b>               | <b>A-dec 300 Continental</b>                                 | <b>A-dec 300 Traditional</b>                                 |
|-------------------------------------|--------------------------------|--|--|
| Control block                       | Standard 4-position            | Standard 3-position, optional 4-position                     | Standard 3-position, optional 4-position                     |
| Left/right conversion               | Yes                            | Yes  | Yes  |
| Tubing options                      | Vinyl or silicone              | Vinyl or silicone  | Vinyl or silicone  |
| Handpiece holder positions          | 4                              | 4  | 4  |
| Brake handle(s)                     | N/A                            | 1 standard, 2nd optional                                     | 1 standard, 2nd optional                                     |
| Touchpad options                    | None, standard or deluxe       | None, standard or deluxe                                     | None, standard or deluxe                                     |
| Quad-voltage intraoral light source | Optional                       | Optional   | Optional   |
| Foot control                        | Disc with chip blower or lever | Disc, disc with chip blower or lever                         | Disc, disc with chip blower or lever                         |
| Arm system                          | Balanced flexarm               | Balanced flexarm with air brake<br>A-dec 335 - support mount | Balanced flexarm with air brake<br>A-dec 334 - support mount |





# Building Relationships.

A-dec is headquartered in Newberg, Oregon USA, but our customers span six continents — every one except Antarctica. Schools and institutions in more than 85 countries rely on A-dec equipment to withstand the rigors of continual use day after day, year after year. Yet with A-dec, you get more than quality equipment. You get full service and support long after the sale and installation. You get local dealers and service professionals, who will respond to your needs and keep you working. And, you get the benefit of more than 50 years of experience to help your practice grow. If you ever have a question, concern or even an idea for improvement, we're just a phone call away. We're here to help you succeed, because at A-dec we don't just build equipment. We build relationships.



#### **A-dec Headquarters**

2601 Crestview Drive  
Newberg, Oregon 97132 USA  
Tel: 1.800.547.1883 within USA/CAN  
Tel: +1.503.538.7478 outside USA/CAN  
Fax: +1.503.538.0276  
a-dec.com

#### **A-dec Australia**

Tel: 1.800.225.010 within AUS  
Tel: +61.(0)2.8332.4000 outside AUS

#### **A-dec China**

Tel: 400.600.5434 within CN  
Tel: +86.571.89026088 outside CN

#### **A-dec United Kingdom**

Tel: 0800.233285 within UK  
Tel: +44.(0).24.7635.0901 outside UK

A-dec, the A-dec logo, A-dec Inspire, Cascade, Century Plus, Continental, Decade, ICX, ICV, Performer, Preference, Preference Collection, Preference ICC, Radius, and reliablecreativesolutions are trademarks of A-dec Inc. and are registered in the United States and other countries. A-dec 500, A-dec 400, A-dec 300, and A-dec 200 are also trademarks of A-dec Inc. None of the trademarks or trade names in this document may be reproduced, copied, or manipulated in any manner without the express, written approval of the trademark owner.

Certain touchpad symbols are proprietary to A-dec Inc. Any use of these symbols, in whole or in part, without the express written consent of A-dec Inc. is strictly prohibited.

The colors and products shown may vary from the actual materials and are subject to change without notice. Please see your authorized A-dec dealer for accurate color samples and the most current product information.

©2017 A-dec Inc. All rights reserved.  
85.0093.00/OP/3M/8-17(11-17)/Rev D2

